SIMPLERA™

APP USER GUIDE

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Medtronic

Medtronic and Medtronic logo are trademarks of Medtronic. ^{TM*} Third-party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

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Medtronic Australasia Pty. Ltd.
Tel: 1800 668 670

Bangladesh:

Sonargaon Healthcare Pvt Ltd.
Mobile: (+91)-9903995417
or (+880)-1714217131

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Medtronic Medikal Teknoloji
Ticaret Ltd. Sirketi.
Tel: +90 216 4694330

USA:

Medtronic Diabetes Global
Headquarters
24-Hour Technical
Support: +1-800-646-4633
To order supplies: +1-800-843-6687

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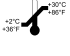
Simplera™ app

Introduction

The Simplera app helps patients manage diabetes. The Simplera app has many key capabilities:

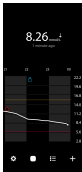
- It records glucose values throughout the day and night.
- It displays glucose values and sends alerts for glucose events in a convenient and discreet manner on a mobile device.
- It displays a history of diet, exercise, and medication events logged by the patient.
- It is compatible with the Apple Watch™* which can display glucose values and receive status and sensor glucose (SG) alerts.

Note: The Simplera app should only be used with supported mobile devices and operating systems. Refer to the local Medtronic website or the local Medtronic support representative for information about compatible mobile devices and operating systems. The Simplera app may not be available in all regions.

Note: This user guide contains some instructions that are specific to the Android™** or iOS™** platform. Where applicable, follow the instructions specific to the platform in use.

System description

The Simplera system includes the following components:

Simplera system	
	Simplera app for iOS™** (MMT-8400) or Simplera app for Android™** (MMT-8401)
	Simplera sensor (MMT-5100J)

Continuous glucose monitoring (CGM) is a technology that lets patients continuously view glucose values. The Simplera app uses an inserted Simplera sensor with Bluetooth™** wireless technology that converts small amounts of glucose from the interstitial fluid under the skin into an electronic signal. The Simplera sensor uses that signal to provide SG values to the Simplera app on a compatible mobile device.

The Simplera app also provides alerts based on SG levels and contains a user interface for entering data such as insulin, meals, exercise, blood glucose (BG) values, and uploading information to the CareLink Personal website.

The Simplera app is available to download from the Apple™ App Store™ or the Google Play™ store and requires the Simplera sensor to function.

Intended use

The Simplera app in combination with the Simplera sensor is intended for use only by patients and caregivers using a compatible mobile device and operating system, and who have sufficient experience to adjust mobile device audio and notification settings. The system is intended to communicate with digitally connected devices for the purpose of managing diabetes.

The Simplera system is designed to replace fingerstick BG testing for diabetes treatment decisions. Treatment decisions should be made based on the combination of SG readings and trend arrows. The system can track trends in glucose concentrations, and aid in detecting events of hypoglycemia and hyperglycemia and helping both current and long-term therapy decisions.

The Simplera system is intended for use in home and professional healthcare environments.

Indications for use

The Simplera system is a real-time CGM system indicated for the management of diabetes in persons ages 2 years and older.

Intended target population

The intended target population for the Simplera system includes children ages 2 years and older, adolescents, and adults.

Intended users

The Simplera system is intended for personal use by individuals to assist in the management of their diabetes, or for use by parents and caregivers who assist these individuals with diabetes management.

Contraindications

The Simplera system has no known contraindications.

Intended clinical benefits

The Simplera system provides information that is used for diabetes management but does not provide any direct therapy. Specifically, the low and high glucose alerts provided by the app when used in combination with the Simplera sensor may allow the user to take appropriate actions to prevent or minimize the severity of hypoglycemia and hyperglycemia.

User safety

Warnings

General

- CGM is not recommended for people who are unwilling or unable to perform blood glucose tests as required or for people who are unwilling or unable to maintain contact with their healthcare professional.
- Always consult a healthcare professional before using sensor glucose values to make treatment decisions if a medication that contains acetaminophen or paracetamol is taken while wearing the sensor. Medications that contain acetaminophen or paracetamol can falsely raise sensor glucose readings. The level of inaccuracy depends on the amount of

acetaminophen active in the body and can differ for each person. Falsely elevated sensor readings can result in over-administration of insulin, which can cause hypoglycemia.

- Medications that contain acetaminophen or paracetamol include, but are not limited to, cold medicines and fever reducers. Check the label of any medications being taken to see if acetaminophen or paracetamol is an active ingredient. Use additional blood glucose meter readings to confirm blood glucose levels.
- Do not use continuous glucose monitoring if hydroxyurea, also known as hydroxycarbamide, is taken. Hydroxyurea is used to treat certain diseases, such as cancer and sickle cell anemia. Hydroxyurea use results in higher sensor glucose readings compared to blood glucose readings. Taking hydroxyurea while using continuous glucose monitoring can result in inaccurate or missed alerts, and substantially higher sensor glucose readings in reports than actual blood glucose readings.
- Always check the label of any medication being taken to confirm if hydroxyurea or hydroxycarbamide is an active ingredient. If hydroxyurea is taken, consult a healthcare professional. Use additional blood glucose meter readings to verify glucose levels.
- If a serious incident related to the device occurs, immediately report the incident to Medtronic and to the applicable competent authority having jurisdiction in their locale.

App and mobile device

- Follow the instructions and safety warnings in this user guide to receive alerts. Missing alerts from the Simplera app can result in undetected low and high glucose levels.
- Always allow notifications for the Simplera app. If notifications are turned off, the app will not send any alerts, including Urgent Low Alert.
- Do not use the Simplera app without understanding how the mobile device settings work. If the mobile device settings are not set up correctly, the app may not send sensor glucose alerts.
- Do not use the Simplera app without sufficient vision and hearing, to allow recognition of the alerts sent by the app.
- Make sure Bluetooth is on, even if the mobile device is in airplane mode. If Bluetooth is off, the app will not send sensor glucose information or alerts.
- Do not use the Simplera app if the mobile device screen or speakers are damaged. If the mobile device is damaged, the app may not send sensor glucose alerts and sensor glucose information may not be shown correctly.
- Do not force close the Simplera app. If the app is closed, the app will not send sensor glucose information or alerts.
- Check the Simplera app occasionally to make sure it is running. The mobile device may close the Simplera app automatically when another app is in use, such as a game. If the Simplera app is closed, the app will not send sensor glucose alerts.
- Do not let the mobile device shut down due to low battery, or the app will not send sensor glucose alerts. Use of the app may deplete the mobile device battery more quickly. Have a charger available to charge the battery if needed.
- Always make sure to open the app after the mobile device restarts to ensure the app sends sensor glucose alerts.
- Make sure to set the snooze to a short enough time so that the app sends an alert again if glucose levels do not improve. When a sensor glucose alert is snoozed, the app will not send that alert again during the length of the set snooze time.
- Do not root or jailbreak the mobile device. Rooting the Android device or jailbreaking the iOS device means to change the software in a way the manufacturer did not intend. If the mobile device is changed in this way, the Simplera app will display an error message when launched and will not continue to operate.

- Turning off automatic updates on the mobile device may help to avoid unintentionally updating to an operating system that is not confirmed as compatible with the Simplera app.
- Always protect the mobile device with a passcode, face authentication, or fingerprint authentication. Sharing the passcode could compromise the security of the device.
- When power saving settings are enabled, the mobile device may delay the alerts and notifications from the Simplera app.

Android users

- Allow **Do Not Disturb Permissions** and **Notifications** for the Simplera app. If **Do Not Disturb Permissions** or **Notifications** are turned off, the app will not send any alerts, including Urgent Low Alert.
- The **Digital Wellbeing** feature is intended to reduce alerts. If this feature is enabled, the Simplera app will not send alerts. If the app timer is set in the **Digital Wellbeing** feature for the Simplera app, the app will shut down and will not send any sensor glucose alerts when the timer expires.
- Do not lower the vibration level to the lowest setting in the Android **Vibration intensity** menu. If vibration is lowered to the lowest setting, the app will not send any vibrations with the alerts, including Urgent Low Alert.
- Alerts for the Simplera app will sound through headphones when headphones are connected. If headphones are connected but they are not being used, sensor glucose alerts may not be heard.

iOS users

- Allow **Critical Alerts** and **Notifications** for the Simplera app. If **Critical Alerts** or **Notifications** are turned off, the app will not send any alerts, including Urgent Low Alert.
- The **Downtime** setting within the **Screen Time** feature is intended to prevent alerts during the **Downtime** period. If this setting is enabled, the Simplera app will be prevented from sending alerts.
- Do not turn off vibration in the iOS Accessibility menu. If vibration is turned off, the app will not send any vibrations with the alerts, including Urgent Low Alert. Vibrate on silent must be turned on in the iOS phone settings or the alerts, including Urgent Low, will not vibrate during Mute All Alerts.

Apple Watch Users

- Always use a blood glucose meter reading to make treatment decisions if no sensor data are available or if symptoms do not match the sensor glucose value.
- To receive sensor glucose information or alerts on the Apple Watch, ensure Bluetooth is enabled and the watch is within range of the mobile device.
- If the mobile device is damaged, the watch may not receive sensor glucose alerts. Do not use the watch if the screen or speakers are damaged.
- The Apple Watch will not vibrate for alerts from the Simplera app if **Haptic Alerts** are disabled.
- Do not disable the mirror feature for the Simplera app in the Apple Watch settings. If the mirror feature is disabled, the watch will not display any alerts, including Urgent Low Alert.
- If the mobile device is unlocked, sensor glucose alerts are not sent to the Apple Watch. Refer to the Simplera app on the mobile device for sensor glucose alerts.
- The Apple Watch may display connection requests. Do not accept connection requests from unknown devices.
- Always protect the Apple Watch with a passcode. Sharing the passcode could compromise the security of the device.

Sensor

- Read this entire user guide before attempting to insert the Simplera sensor. The inserter portion of the sensor does not work the same way as other Medtronic insertion devices. The sensor is not inserted the same way as other Medtronic sensors. Failure to follow directions may result in improper insertion, pain, or injury.
- Do not use continuous glucose monitoring if hydroxyurea, also known as hydroxycarbamide, is taken. Hydroxyurea is used to treat certain diseases, such as cancer and sickle cell anemia. Hydroxyurea use results in higher sensor glucose readings compared to blood glucose readings. Taking hydroxyurea while using continuous glucose monitoring can result in substantially higher sensor glucose readings in reports than actual blood glucose readings.
- Always check the label of any medication being taken to confirm if hydroxyurea or hydroxycarbamide is an active ingredient. If hydroxyurea is taken, consult a healthcare professional. Use additional blood glucose meter readings to verify glucose levels.
- Taking medications that contain acetaminophen or paracetamol, including, but not limited to fever reducers and cold medicine, while wearing the sensor, may falsely raise SG readings. The level of inaccuracy depends on the amount of acetaminophen or paracetamol active in the body and may be different for each person. Always check the label of any medications to confirm whether acetaminophen or paracetamol is an active ingredient.
- Always examine the Simplera sensor box for damage. If the sensor box is open or damaged, examine the sensor for damage. If the sensor is visibly damaged, discard the device to avoid possible contamination.
- Do not use the Simplera sensor if any part of the device is damaged. If the device is damaged, discard the device to avoid possible contamination.
- Do not use the Simplera sensor if the tamper band is broken, damaged, or is missing from the device. The sensor is sterile and non-pyrogenic unless the device is damaged. If the tamper band is broken, damaged, or is missing from the device, the sensor and needle can be exposed to contamination. A sensor and needle exposed to contamination can cause insertion site infection if inserted into the body.
- Do not use the Simplera sensor if the cap label is broken, damaged or missing from the device. The sensor is sterile and non-pyrogenic unless the device is damaged. If the cap label is broken, damaged, or missing from the device, the sensor and the needle can be exposed to contamination. A sensor and needle exposed to contamination can cause insertion site infection if inserted into the body.
- Do not unscrew or remove the Simplera sensor cap until the device is ready to be used. Do not remove the cap and store the device for future use. The sensor is sterile and non-pyrogenic unless the cap is removed from the device or the tamper band is broken. If the cap is not on the device or the tamper band is broken, the sensor and the needle can be exposed to contamination. A sensor and needle exposed to contamination can cause insertion site infection if inserted into the body.
- Do not remove the cap and place it back on the device. Placing the cap back on the device could cause damage to the needle, prevent a successful insertion, and cause injury.
- Do not change or modify the Simplera sensor. Changing or modifying the sensor can result in improper insertion, pain, or injury.
- Do not let children hold the Simplera sensor without adult supervision. Do not let children put any part of the Simplera sensor in their mouth. This product poses a choking hazard for young children that can result in serious injury or death.
- Watch for bleeding at the insertion site on top of the Simplera sensor. If bleeding occurs, apply steady pressure with a sterile gauze pad or clean cloth placed on top of the sensor for up to three minutes. If bleeding continues, is significantly visible on top of the sensor, or if there is excessive pain or discomfort after insertion, follow these steps:

1. Remove the Simplera sensor and continue to apply steady pressure until the bleeding stops.
 2. Dispose of the Simplera sensor. See *Disposal*, page 47.
 3. Check the site for redness, bleeding, irritation, pain, tenderness, or inflammation. If there is redness, bleeding, irritation, pain, tenderness, or inflammation, contact a healthcare professional.
 4. Insert a new Simplera sensor in a different location.
- Adverse reactions can cause serious injury.

Precautions

Diabetes treatment decisions should be made based on a combination of SG readings, trend arrows, glucose target ranges, active alerts, and recent events (such as insulin doses, exercise, meals, and medications).

Using CGM information to make treatment decisions

After becoming familiar with CGM, treatment decisions should be made based on all the information available, including the following:

- SG readings
- Trend arrows
- Active SG alerts
- Recent events such as insulin doses, medication, meals, exercise, etc.

Consult a healthcare professional to determine the right glucose target ranges.

Sensor

- Do not use the Simplera sensor adjacent to other electrical equipment that may cause interference with normal system operation. For more information on electrical equipment that may cause interference with normal system operation, see *Exposure to magnetic fields and radiation*, page 18.
- Some skin care products, such as sunscreens and insect repellents, can damage the Simplera sensor. Do not allow skin care products to touch the sensor. Wash hands after using skin care products before touching the sensor. If any skin care products touch the sensor, immediately wipe the sensor with a clean cloth.

Exposure to magnetic fields and radiation

- Do not expose the Simplera sensor to Magnetic Resonance Imaging (MRI) equipment, diathermy devices, or other devices that generate strong magnetic fields (for example, CT scan, or other types of radiation). Exposure to strong magnetic fields can cause the sensor to malfunction, result in serious injury, or be unsafe.

IEC 60601-1-2; Special EMC Precautions for Medical Electrical Equipment

1. Special Precautions regarding Electromagnetic Compatibility (EMC): This body worn device is intended to be operated within a reasonable residential, domestic, public or work environment where common levels of radiated "E" (V/m) or "H" fields (A/m) exist, such as cellular phones, Wi-Fi™*, Bluetooth® wireless technology, electric can openers, microwave and induction ovens. This device generates, uses, and can radiate RF energy and, if not installed and used in accordance with the provided instructions, may cause harmful interference to radio communications.

2. Portable and mobile RF communications equipment can affect medical electrical equipment. If you encounter RF interference from a mobile or stationary RF transmitter, move away from the RF transmitter that causes the interference.
3. Be careful when using the sensor closer than 30 cm (12 in) to portable RF equipment or electrical equipment. If the sensor must be used next to portable RF equipment or electrical equipment, observe the sensor to verify correct system operation. Degradation of the performance of the sensor could result.

Risks and side effects

General risks with the sensor include the following:

- Skin irritation or other reactions
- Bruising
- Discomfort
- Redness
- Bleeding
- Pain
- Rash
- Infection
- Raised bump
- Appearance of a small “freckle-like” dot where the sensor was inserted
- Allergic reaction
- Fainting secondary to anxiety or fear of needle insertion
- Soreness or tenderness
- Swelling at insertion site
- Sensor filament breakage or damage
- Minimal blood splatter associated with sensor needle removal
- Residual redness associated with adhesive or tapes or both
- Scarring

Failure to use the Simplera app according to the instructions for use may result in hypoglycemia or hyperglycemia. If your glucose alerts and sensor glucose readings do not match your symptoms, or if no sensor data are available, use a fingerstick blood glucose value from your blood glucose meter to make diabetes treatment decisions. Seek medical attention when appropriate. Consult a healthcare professional about how to use the information displayed on the app to help manage your diabetes.

Hazardous substances

For materials information such as compliance with European Union (EU) Registration, Evaluation, Authorization and Restriction of Chemicals (REACH) regulation, EU Restriction of Hazardous Substances (RoHS) directive, and other product stewardship program requirements, please visit <http://www.medtronic.com/productstewardship>.

Allergens

The Simplera sensor contains nickel in stainless steel.

Reagents

The sensor contains two biological reagents: glucose oxidase and human serum albumin (HSA). Glucose oxidase is derived from *Aspergillus niger* and manufactured to meet industry

requirements for the extraction and purification of enzymes for use in diagnostic, immunodiagnostic, and biotechnical applications. The HSA used on the sensor consists of purified and dried albumin fraction V, derived from pasteurized human serum, which is cross-linked via glutaraldehyde. Approximately 3 µg of glucose oxidase and approximately 10 µg of HSA are used to manufacture each sensor.

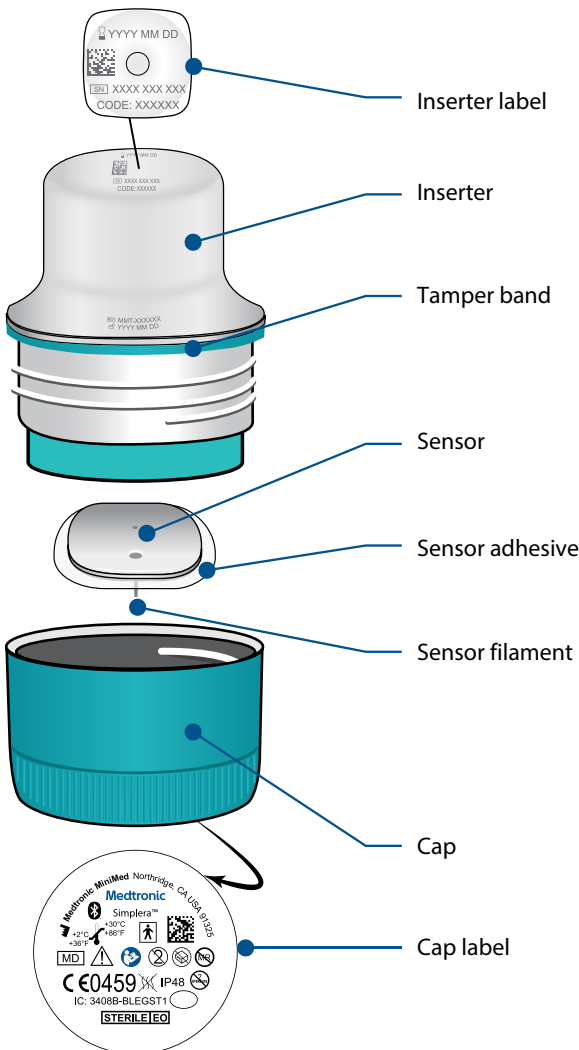
Assistance

Contact the local representative or refer to the local Medtronic website for assistance. Refer to the Medtronic Diabetes International Contacts list in this user guide for contact information.

Simplera app setup

Search for and download the Simplera app  from the Apple App Store or Google Play store on the supported mobile device. To set up the Simplera app, follow the instructions on the screen.

Simplera sensor device components


New sensor setup

Follow the video instructions on how to insert the sensor. Then follow the on-screen instructions to complete the sensor start up.


Where to insert the Simplera sensor

The images that follow show insertion sites for ages 2 to 17 years, and ages 18 years and older. Choose an insertion site for the applicable age group. Target the shaded areas shown in the image, and make sure that the insertion site has a sufficient amount of fat.

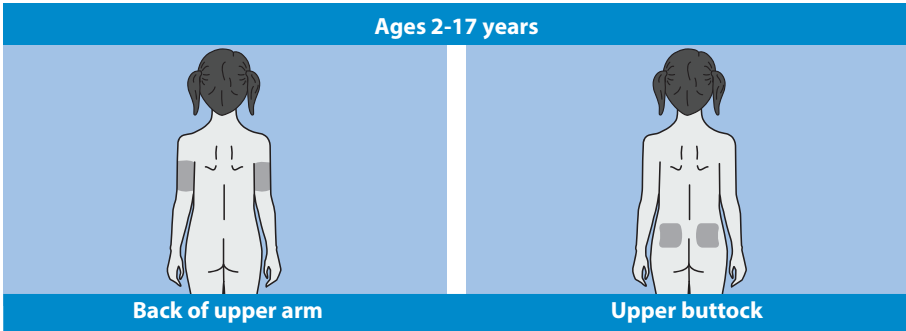
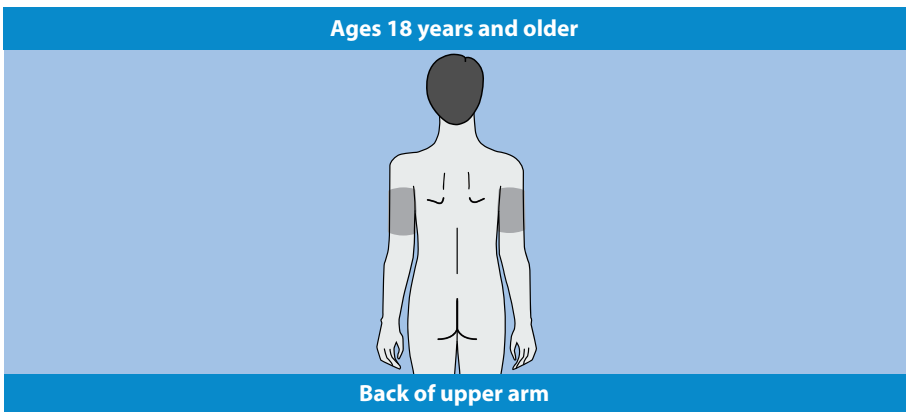
Insertion on abdomen for ages 2 to 17 years has not been evaluated for accuracy.


Note: Insertion into the upper buttocks should target the top-third of the buttocks area. Assistance from another person may be needed for sensor insertion into the back of the upper arm or the upper buttock. If assistance is not needed, a mirror may be helpful for self-insertion.

Insertion on abdomen and upper buttocks for ages 18 years and older has not been evaluated for accuracy.

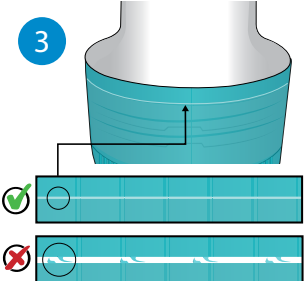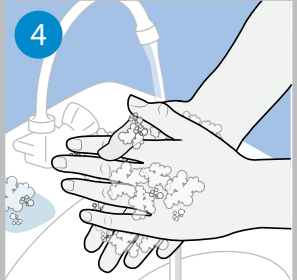
Inserting the Simplera sensor

Preparing for insertion	
<p>1</p> 	<p>The inserter label is on the top of the inserter.</p> <p>1. Before insertion, perform the following steps:</p> <ul style="list-style-type: none">• Inspect the expiration date. Do not use an expired Simplera sensor.• Make note of the serial number (SN) and the CODE. Both numbers will be used later to pair the sensor with the Simplera app. <p>Note: The SN and CODE label is also on the inside of the Simplera sensor box lid.</p>
<p>2</p> 	<p>2. Inspect the cap label for damage before insertion.</p> <p>Note: Do not use the Simplera sensor if the cap label is damaged or missing from the cap.</p>
<p>3</p> 	<p>3. Inspect the tamper band to make sure that it is not broken, damaged, or missing from the device.</p> <p>Note: Do not use the Simplera sensor if the tamper band is broken, damaged or missing.</p>
<p>4</p> 	<p>4. Wash hands thoroughly with soap and water.</p> <p>Note: Wear gloves when inserting the Simplera sensor into another person to avoid accidental contact with patient blood. Minimal bleeding can occur.</p>

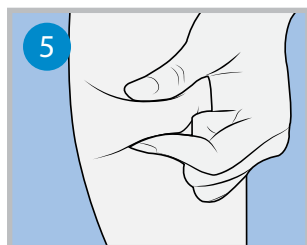
5. Choose an insertion site that has a sufficient amount of fat. For insertion sites see *Where to insert the Simplera sensor*, page 21. For the best sensor performance, and to avoid accidental sensor removal, do not insert the Simplera sensor into the following areas:

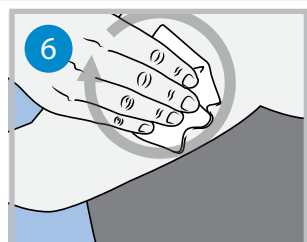- muscle, tough skin, or scar tissue
- areas that are constrained by clothing or accessories
- areas subjected to rigorous movement during exercise

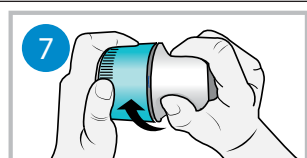
6. Clean the insertion site with alcohol. Allow the insertion site to air dry.

7. Unscrew the cap from the inserter, breaking the tamper band.

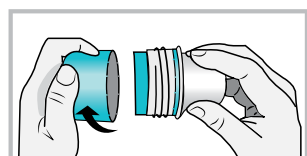
Note: Do not use the Simplera sensor if the tamper band is broken, damaged, or missing.


Insertion

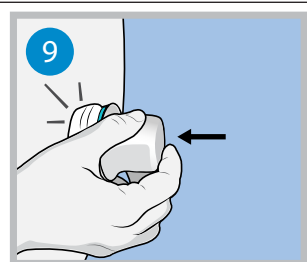8. Place the inserter on top of the prepared insertion site.

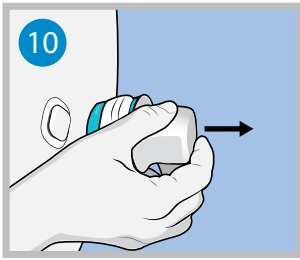
9. Press the inserter firmly against the body until there is a click.


10. Gently pull the inserter straight from the body.

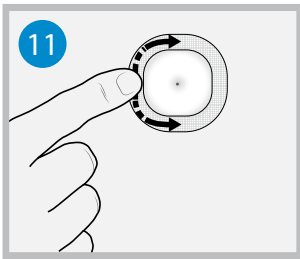After insertion

11. Smooth down the sensor adhesive with a finger to ensure the sensor stays on the body for the entire duration of use.

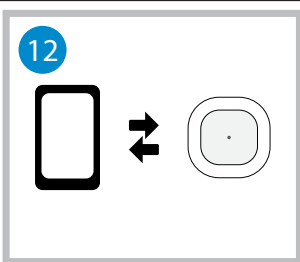Note: Use over-the-counter tape if desired for additional adhesion.

12. Pair the Simplera sensor with a compatible mobile device.

Note: The SN and CODE are required to pair the sensor with a compatible mobile device. For details on how to pair the sensor with a compatible mobile device and the Simplera app, see *Pairing the Simplera sensor, page 24*.

Pairing the Simplera sensor

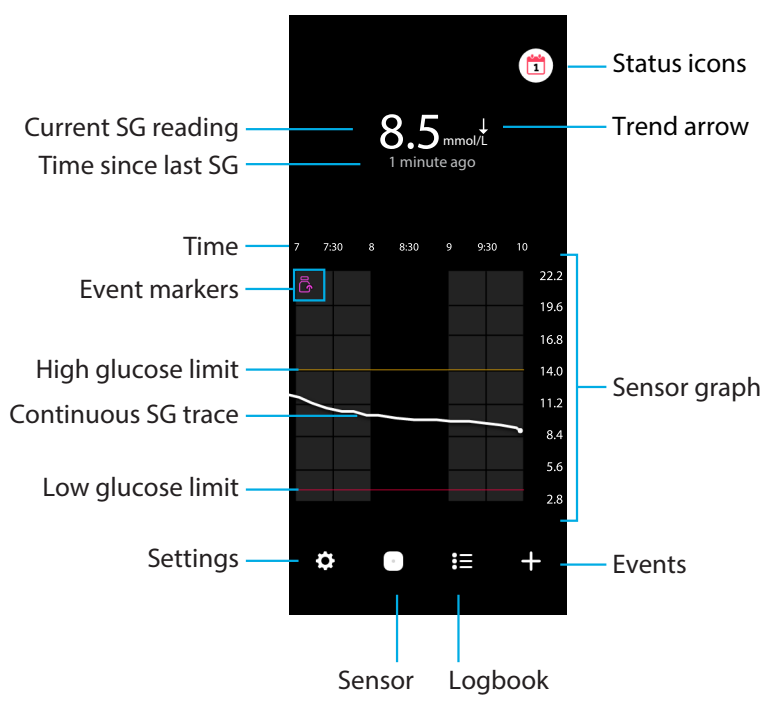
The Simplera sensor must be paired to the compatible mobile device before use. Follow the on-screen instructions to pair the sensor to the mobile device.

Completing the app setup

Continue to follow the on-screen instructions to enable notifications. For more information, see *Glucose alerts information, page 32*.

Home screen

The following figure shows the Home screen of the app.

Note: This screen may vary depending on the mobile device and platform.


Note: The Home screen does not display sensor information the first time the app is opened. The first SG reading appears after the sensor is successfully paired and warm-up is complete.

Item	Description
Current SG reading	Displays the current SG reading, which the sensor calculates and sends wirelessly to the app.
Trend arrows	Display the glucose trend and the rate at which the most recent SG level has risen or fallen. For details on how to set the rate of change for falling and rising SG values, see <i>Glucose alerts</i> , page 33.
Time since last SG	Displays the time since the last SG reading that the sensor has sent wirelessly to the app.
Event marker	Displays an icon for specific events, such as exercise, BG readings entered, carbs eaten, or insulin taken.
Continuous SG trace	Displays current and previous SG readings. Tap any point on the graph to view the details of the selected SG reading or event. For more navigation tips, see <i>Sensor graph</i> , page 37.
Status icons	Display the corresponding status icon if the sensor life has expired, the connection with the sensor is lost, or alerts are muted. For details, see <i>Sensor screen</i> , page 28.
Time	Displays preset time intervals of 3, 6, 12, and 24 hours.

Item	Description
High glucose limit and low glucose limit	Displays a line indicating the high and low SG limits on the sensor graph. The orange line indicates the high SG limit; the red line indicates the low SG limit.
Sensor graph	Displays a white line representing SG readings over a selected time interval. It also indicates high and low SG limits. For details on the sensor graph, see <i>Sensor graph, page 37</i> .
Settings	Displays the settings screen to customize additional features, such as alert volume, mute all alerts, glucose alerts, CareLink, units, and provides access to the user guide, getting started guide, and help.
Sensor	Displays the Simplera sensor status information and allows the user to replace a sensor or unpair a sensor. For details, see <i>Sensor screen, page 28</i> .
Logbook	Displays up to 30 days of events, and alerts in chronological order.
Events	Displays the events screen to enter additional information for exercise, carbs, or insulin taken. For details, see <i>Sensor screen, page 28</i> .

Settings screen

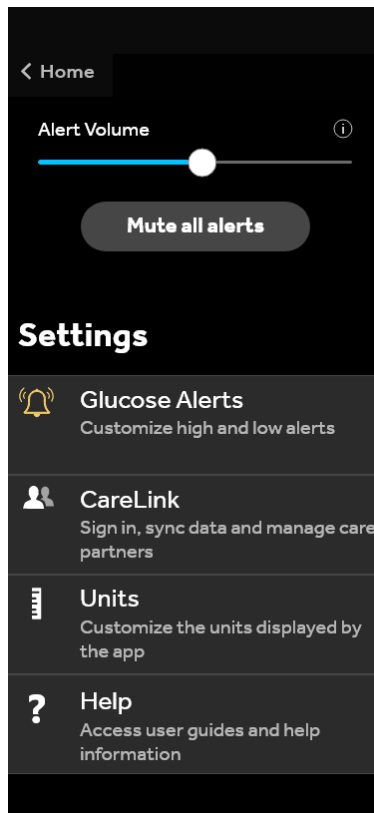
Tap  to see the Settings screen.

Item	Description
Home	Closes the Settings screen and returns to the Home screen.
Alert Volume	Allows the option to adjust the volume of alerts. Alerts will sound louder and repeat if not acknowledged. For details, see <i>Glucose alerts information, page 32</i> .
Alert Volume information	Displays a description of Alert Volume slider and Mute All Alerts button. For details, see <i>Glucose alerts information, page 32</i> .
Mute all alerts	Allows the option to mute all alerts for a specified duration. Urgent Low Alert will always vibrate. For details, see <i>Glucose alerts information, page 32</i> . Note: Urgent Low Alert will not vibrate if vibration is turned off in the iOS Accessibility menu or lowered to the lowest setting in the Android Vibration intensity menu.
Glucose Alerts	Display the options to set glucose monitoring alerts. For details, see <i>Glucose alerts information, page 32</i> .
CareLink	Provides access to a menu to change the CareLink user login information. For details, see <i>Glucose alerts information, page 32</i> .
Units	Allows the option to display glucose in mg/dL or mmol/L. Note: This feature may not be available in all regions.
Help	Provides the user with software version information, access to the user guide, getting started guide, setup wizard, and the end user license agreement.

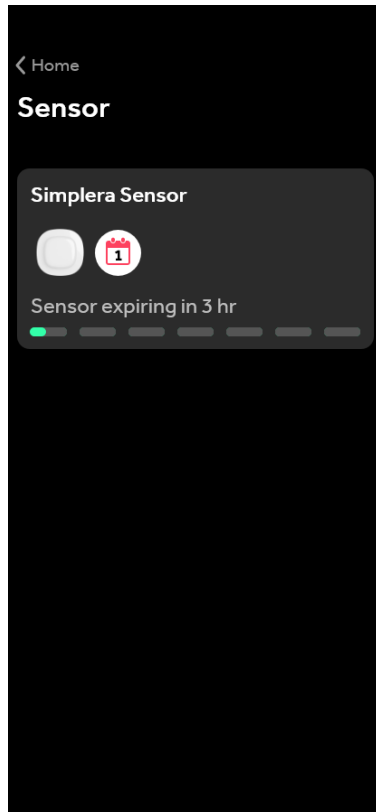Status icons


The status icons appear on the Sensor screen. These icons provide a way to quickly check the status of the sensor. If any condition becomes critical and requires immediate attention, the corresponding status icon is also displayed on the Simplera app Home screen. The icons displayed on the Home screen are interactive and provide more information about the current status.


Icon	Icon name	Description
	Sensor life	Indicates there is fewer than one day left of sensor life.
	Sensor expired	The sensor has expired and must be replaced.
	Communication error	Indicates there is a communication error, the sensor is not paired to the mobile device, or Bluetooth™ is off.
	General error	Indicates a general error.

Sensor screen

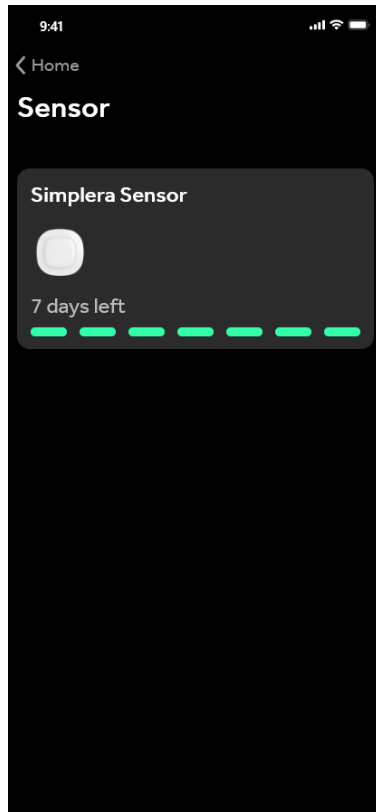
Tap  at the bottom of the Home screen to see the Sensor screen. This screen displays the sensor status as well as status icons.

Tap  to see sensor information.

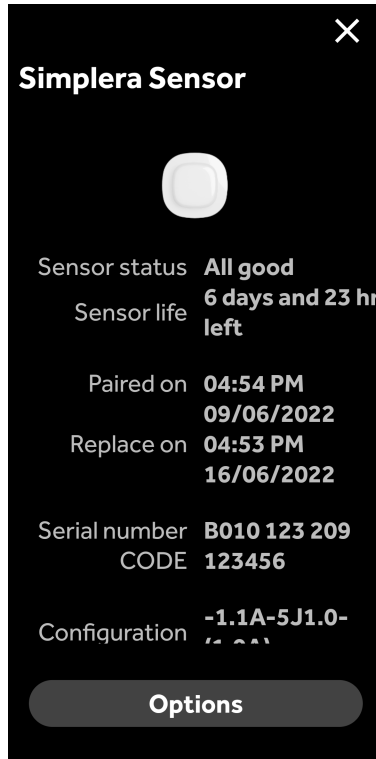Tap **Options** to show the menu.

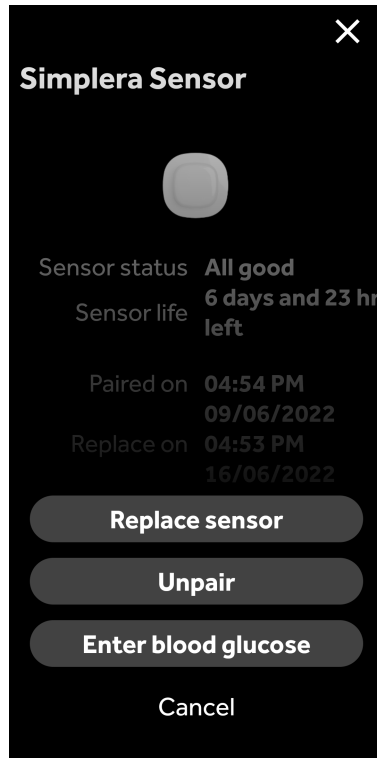Item	Description
Sensor status	Sensor events or error messages are displayed.
Sensor life	The time remaining until the sensor must be replaced is displayed.
Paired on	The time and date the sensor was last paired is displayed.
Replace on	The estimated time and date when the sensor must be replaced is displayed.
Replace sensor	Allows the addition of a new sensor.
Unpair	Provides directions to unpair a sensor from the mobile device settings.
Enter blood glucose	Shows the screen that allows entry of blood glucose values.
Cancel	Returns to the previous screen.

Note: The Bluetooth devices list in the mobile device will grow as new sensors are paired. Periodically review the list and remove old sensors. Avoid removing the currently paired sensor from the list, or it will need to be paired again in the app. The current sensor will appear in the paired devices list, or it will show that it is connected.

Enter blood glucose as an optional calibration

While the Simplera system does not require calibration, the Simplera app gives the option to enter BG meter readings as an optional calibration. Enter the BG meter reading into the app immediately after checking BG. Review the BG meter instructions for guidance on how to check the BG.

- Always make sure hands are clean before testing BG.
- Avoid use of an old BG reading or reuse of BG readings from previous calibrations.
- If for some reason calibration is unsuccessful, wait at least 15 minutes before attempting another calibration.

To enter an optional calibration into the Sensor screen:

1. Measure the BG with a BG meter.
2. Tap  in the lower center of the Home screen to see the Sensor screen.
3. Tap on the sensor card to see the sensor details.
4. Tap **Options** to show the menu.
5. Tap **Enter blood glucose**.
6. Use the number pad to enter the BG meter reading.
7. Make sure the value displayed above the number pad is correct. If the value is incorrect, tap  to clear it and enter the correct value.
8. To calibrate the sensor with the entered BG meter reading, tap **Calibrate**. The app returns to the Home screen, and  appears on the graph at the selected time.

Alerts

The Simplera app provides **glucose alerts** and **status alerts**. These alerts provide information about glucose levels and the status of the Simplera system.

Glucose alerts information

Set **glucose alerts** for the following situations:

- SG values are rising or falling faster than the set limits
- SG values have gone above or below the set limits
- SG values are predicted to go above or below the set limits

The **glucose alerts** in the Simplera app are listed in the following table:

Alert type	Description
High Alert	SG values have gone above the set high limit.
High predicted	The SG value is predicted to go above the set high limit, within a set period of time (up to 60 minutes ahead).
Rise Alert	The SG value is rising faster than the set rate (corresponding to the rising arrows displayed next to the SG level).
Low Alert	The SG value has gone below the set low limit.
Low predicted	The SG value is predicted to go below the set low limit, within a set period of time (up to 60 minutes ahead).
Fall Alert	The SG value is falling faster than the set rate (corresponding to the falling arrows displayed next to the SG level).
Urgent Low Alert	The SG level is at or below 3.0 mmol/L (54 mg/dL).

Note: If the Simplera app detects that alerts are disabled, the app will not send any alerts, including the Urgent Low Alert.

Glucose alerts can be customized and are set up as described in *Glucose alerts*, page 33.

Status alerts

The Simplera app also has **status alerts** that provide information about actions needed to ensure the correct functioning of the system. See the status alerts table below for a complete listing of these alerts. For more information on how to address these alerts, see *Troubleshooting*, page 44.

The **status alerts** in the Simplera app are listed below:

Make sure to act on these alerts to continue to receive SG information in the future.

Alert type	Description
Mobile device battery low	The mobile device battery has reached or fallen below 20% of its power.
Sensor warmup	The sensor is warming up. This will take 2 hours.
Calibration not accepted	The BG meter value could not be used to calibrate.
Change sensor	The sensor is not working properly and needs to be replaced.
Lost communication	The app and sensor have not been communicating for 30 minutes. The app may have closed if there are too many apps running at the same time or if there is radio frequency interference.
Insufficient battery	The sensor has insufficient battery. Change the sensor.
Sensor end of life	The sensor has reached its maximum life of 7 days. The Simplera app is no longer receiving sensor information. Insert a new sensor.
Sensor updating	The sensor is updating. Updating can take up to 2 hours. SG readings will not be available during this time. Measure the BG value with the meter.
Connection error	The sensor is trying to reconnect. Wait at least 30 minutes.
Jailbroken device detected	The software on the iOS mobile device has been changed such that it no longer works in the way the manufacturer intended. The app cannot be used with a jailbroken device. Changing the mobile device operating software causes the app to stop working.
Rooted device detected	The software on the Android mobile device has been changed such that it no longer works in the way the manufacturer intended. The app cannot be used with a rooted device. Changing the mobile device operating software causes the app to stop working.

Glucose alerts

Setting low alerts

Low alert settings include the following:

Low alert setting	Description
Day Starts at	Set when the day begins, what the low limit is during the day, and when the app should send alerts at this setting.
Night Starts at	Set when the night begins, what the low limit is during the night, and when the app should send alerts at this setting.
Low Limit	The Low Limit is the SG value on which the other low settings are based. The low limit can be set from 3.3 to 5.0 mmol/L (60 to 90 mg/dL). On the SG graph, the low limit appears as a red horizontal line at the set value.
Don't alert me	When Don't alert me is selected, the app will not send any alerts for low glucose. The app still sends an Urgent Low Alert since this alert is always on.

Low alert setting	Description
At Low Limit	When At Low Limit is selected, the app displays a Low Sensor Glucose alert when the SG value reaches or falls below the low limit.
Before Low Limit	When Before Low Limit is selected, the app sends a Low predicted alert any time the SG is predicted to reach the low limit. This alert notifies of potential low glucose levels before they occur.
Before and at Low Limit	When Before and at Low Limit is selected, the app sends a Low predicted alert any time the SG is predicted to reach the low limit and when the SG value reaches or falls below the low limit.
Time Before Low	Time Before Low is available only when Before Low Limit or Before and at Low Limit is selected. Time Before Low determines when the app will send a Low predicted alert, if the SG values continue to decrease at the current rate of change. The option is available to set a time between ten minutes and one hour.
Max Volume at Night	Max Volume at Night is available only in the Night Time section and when Alert Me is selected. Toggle Max Volume at Night on for low alerts to make all low alerts sound at max volume at night.

To set low alerts:

1. On the Home screen, tap  and select **Glucose Alerts > Low Alerts**.
2. Tap **Day starts at** and set the desired start time. Tap **Save**.
3. Tap **Low Limit** and set the desired limit between 3.3 to 5.0 mmol/L (60 and 90 mg/dL). Tap **Save**.
4. To receive alerts when the SG value approaches the low limit, tap **Alert Me** and choose one of the following:
 - a. At Low Limit
 - b. Before Low Limit
 - c. Before and at Low Limit
5. If **Before Low Limit** or **Before and at Low Limit** is selected, the **Time Before Low** appears. Set the desired duration of time to receive a Low predicted alert and tap **Save**.
6. Return to the Low Alerts screen and tap **Night starts at**. Set the desired start time and tap **Save**. Repeat steps 3 through 5 to set night time alerts.
7. Toggle **Max Volume At Night** on or off. When **Max Volume At Night** is toggled on, the low alerts sound at max volume at night.
8. When the alerts are set up, tap  to return to the Glucose Alerts screen.
9. Tap  to return to the settings screen.
10. Tap  to close the menu and return to the Home screen.

Setting high alerts

High alert settings include the following:

High alert setting	Description
Day Starts at	Set when the day begins, what the high limit is during the day, and when the app should send alerts at this setting.
Night Starts at	Set when the night begins, what the high limit is during the night, and when the app should send alerts at this setting.

High alert setting	Description
High Limit	The high limit is the SG value on which the other high settings are based. The high limit can be set from 5.5 to 22.1 mmol/L (100 to 400 mg/dL). On the SG graph, the high limit appears as an orange horizontal line at the set value.
Don't alert me	When Don't alert me is selected, the app does not send any alerts for high glucose.
At High Limit	When At High Limit is selected, the app displays a High Sensor Glucose alert when the SG value reaches or exceeds the high limit.
Before High Limit	When Before High Limit is selected, the app sends a High predicted alert any time the SG value is predicted to reach the high limit. This notifies of potential high glucose levels before they occur.
Before and at High Limit	When Before and at High Limit is selected, the app sends a High predicted alert any time the SG is predicted to reach the high limit and when the SG value reaches or rises above the high limit.
Time Before High	The Time Before High option is available only when Before High Limit or Before and At High Limit is selected. This option determines when the app sends a High predicted alert if the SG values continue to increase at the current rate of change. Set a time between ten minutes and one hour.
Max Volume at Night	Max Volume at Night is available only in the Night time section and when Alert Me is selected. Toggle Max Volume at Night on for high alerts to make all high alerts sound at max volume at night.

To set high alerts:

1. On the Home screen, tap  and select **Glucose Alerts > High Alerts**.
2. Tap **Day starts at** and set the desired start time. Tap **Save**.
3. Tap **High Limit** and set the desired limit between 5.5 to 22.1 mmol/L (100 and 400 mg/dL). Tap **Save**.
4. To receive alerts when the SG value approaches the high limit, tap **Alert Me** and choose one of the following:
 - a. At High Limit
 - b. Before High Limit
 - c. Before and at High Limit
5. If **Before High Limit** or **Before and at High Limit** is selected, the **Time Before High** appears. Set the desired time to receive a High predicted alert and tap **Save**.
6. Return to the High Alerts screen and tap **Night starts at**. Set the desired start time and tap **Save**. Repeat steps 3 through 5 to set night time alerts.
7. Toggle **Max Volume At Night** on or off. When **Max Volume At Night** is toggled on, high alerts sound at max volume at night.
8. When the alerts are set up, tap  to return to Glucose Alerts screen.
9. Tap  to return to the settings screen.
10. Tap  to close the menu and return to the Home screen.

Setting rate alerts

Rate alerts notify when the SG value is rising (Rise Alert) or falling (Fall Alert) equal to or faster than the set rate. These alerts help to understand how glucose levels are affected, for example, by carbs or exercise.

On the Home screen, these fall or rise rates are indicated by arrows, as shown in *Home screen*, page 24. The more arrows, the faster the rate of change.

↑	SG has been rising at a rate of 0.06 mmol/L (1 mg/dL) or more per minute, but less than 0.11 mmol/L (2 mg/dL) per minute.
↓	SG has been falling at a rate of 0.06 mmol/L (1 mg/dL) or more per minute, but less than 0.11 mmol/L (2 mg/dL) per minute.
↑↑	SG has been rising at a rate of 0.11 mmol/L (2 mg/dL) or more per minute, but less than 0.17 mmol/L (3 mg/dL) per minute.
↓↓	SG has been falling at a rate of 0.11 mmol/L (2 mg/dL) or more per minute, but less than 0.17 mmol/L (3 mg/dL) per minute.
↑↑↑	SG has been rising at a rate of 0.17 mmol/L (3 mg/dL) or more per minute.
↓↓↓	SG has been falling at a rate of 0.17 mmol/L (3 mg/dL) or more per minute.

To set rate alerts:

1. On the Home screen, tap  and select **Glucose Alerts > Fall & Rise Alerts**.
2. Toggle **Fall Alert** on.
3. Tap the arrow option with the desired fall rate.
4. Toggle **Rise Alert** on.
5. Tap the arrow option with the desired rise rate.
6. When fall and rise alerts are set up, tap  to return to Glucose Alerts screen.

Setting alert snooze time

The snooze feature gives the option to set a snooze time for alerts. If the alert condition still persists after the set snooze time, the app sends a notification of the alert. The set snooze time for high and rise alerts may be different from the set snooze time for low and fall alerts.

To set the alert snooze time:

1. On the Home screen, tap  and select **Glucose Alerts > Snooze Time**.
2. Tap **Low & Fall Alerts** and set the desired snooze time for low and fall alerts. Tap **Save**.
3. Tap **High & Rise Alerts** and set the desired snooze time for high and rise alerts. Tap **Save**.

Adjusting alert volume

All Simplera app alerts sound at the volume set in the app. Alerts sound louder and repeat if not acknowledged.

To adjust alert volume:

1. On the Home screen, tap .
2. At the top, there is a slider. Swipe left and right to adjust the volume.

Note: If the volume is set at 0%, a popup message appears stating: "Alert volume set at 0%. If you don't respond to the initial alert, it will sound and repeat. Your alerts are NOT muted. To mute alerts, tap **Mute All Alerts**."

Muting alerts

All Simplera app alerts will be muted for the selected time. Urgent Low Alert will always vibrate.

To mute alerts:

1. On the Home screen, tap .
2. Tap **Mute All Alerts** and 4 options appear.
 - a. 30 min
 - b. 1 hour
 - c. 4 hours (max)
 - d. Custom

If Custom is selected, adjust the mute time to the desired length and tap **Mute**.
3. A popup message appears to confirm the alerts are muted. To cancel the mute, tap **Cancel Mute**.

Acting on Simplera app alerts

The Simplera app sends alerts, similar to notifications sent from other apps on the mobile device. All alerts vibrate according to the mobile device vibration settings.

Open the Simplera app to address an alert. Dismissing the notification only removes it from the list of notifications on the mobile device. If the notification is dismissed from the mobile device, but isn't acted on it in the Simplera app, the alert may repeat.

The alert is displayed on the screen when the app is opened. Select **OK** to clear status alerts. For SG alerts, either swipe the alert up, which will snooze the alert for the set snooze period, or swipe the alert down to set a new snooze period.

Sensor graph


The sensor graph displays the current SG reading. It also gives the option to view a history of SG readings and entered events.

Navigating through the SG graph

- Swipe the center of the graph right and left to view historical data and to return to the current graph location.
- Pinch and stretch the center of the graph to zoom in and out on the graph data.
- Tap the graph twice to view the graph at the selected 3-hour, 6-hour, 12-hour, and 24-hour zoom levels.
- Tap the horizontal time axis above the graph to center the selected time point on the graph. This gives the option to view the details of the selected SG reading or event in an information box that appears above the graph. For details, see *Graph information boxes, page 37*.
- Double tap the vertical SG values axis to the right of the graph to return to the current SG reading displayed on the graph.

Graph information boxes

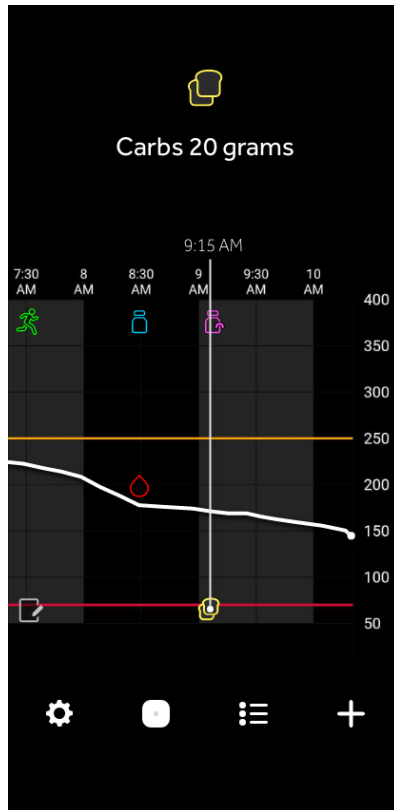Tap on the graph to view information about an event or SG reading. The tapped time is marked by a vertical cursor on the graph, and a box with the information about the tapped event or SG reading appears above the graph. The following figure provides an example:

For SG readings, information boxes display the SG value, its date and time, and trend arrows, similar to the current SG value information displayed on the Home screen. If there is no SG value for the selected time point, the information box displays the sensor status message for that point.

For event markers, information boxes display the details specific to each event type. For details, see *Events screen, page 39*.

Drag the cursor across the graph, to skip at five-minute intervals, to pinpoint a specific SG value on the graph.

Information boxes appear for a few seconds and then close.

Logbook screen

The **Logbook** screen displays a history of alerts and events that occurred on the selected day, with the most recent entries at the top of the list.

To view Logbook entries:

1. On the Home screen, tap  to select **Logbook**.
2. Do any of the following to view the desired information:
 - Tap **alerts** or **events** to filter the list by the specific type. To view the entire list, select **all**.
 - Swipe down and up on the list to view the entries.
 - To delete an event entry, swipe it left and tap **Delete**.

Note: Alerts or calibration events cannot be deleted.

- Tap the desired entry to expand it and view the details.

Events screen

Events help capture information that may affect glucose levels. Use the **Events** screen on the app to enter and save certain types of events.

Event icon	Event name	Description
	Blood Glucose	The Blood Glucose marker represents the BG meter readings.
	Rapid acting insulin	The Rapid acting insulin marker represents the type and amount of insulin delivered.
	Long acting insulin	The Long acting insulin marker represents the type and amount of insulin delivered.
	Carbs	The Carbs marker represents the amount of carbohydrates consumed, either with food or drink.
	Exercise	The Exercise marker represents the intensity and duration of exercise routine.
	Notes	The Notes marker represents other information relevant to diabetes management. For example, a record of other medications taken, illness, or stress.

Note: Make it a practice to enter events when they happen. Events can be deleted and entered again.

Entering BG meter readings

The app gives the option to enter BG meter readings. For example, if measuring BG when eating or when the SG value is rising or falling rapidly, enter the BG meter reading into the app.

Note: Enter a BG value between 0.6 mmol/L and 33.3 mmol/L (10 mg/dL and 600 mg/dL) into the Events screen.

To enter a BG meter reading:

1. Measure the BG with a BG meter.
2. At the bottom of the Home screen, tap .
3. Tap  on the Events screen.

4. Use the number pad to enter the BG meter reading.
5. Make sure the value displayed above the number pad is correct. If the value is incorrect, tap  to clear it and enter the correct value.
6. Tap **Log** at the bottom of the screen. The app returns to the Home screen, and  appears on the graph at the selected time.

To enter an optional calibration, see *Enter blood glucose as an optional calibration, page 31*.

Entering insulin injection information

If delivering insulin using an insulin pump, insulin pen, or a syringe, use the app to record the amount of insulin injected.

To enter the amount of rapid acting insulin injected:

1. Make a note of the amount of insulin taken.
2. At the bottom of the Home screen, tap .
3. Tap  on the Events screen.
4. To change the date or time for the entry, tap **Time** and make the appropriate changes.
5. Use the number pad to enter the insulin amount.
6. Make sure the value displayed above the number pad is correct. If not, tap  to clear it and enter the correct value.
7. Tap **Log** at the bottom of the screen. The app returns to the Home screen, and  appears on the graph at the selected time.

To enter the amount of long acting insulin injected:

1. Make a note of the amount of insulin taken.
2. At the bottom of the Home screen, tap .
3. Tap  on the Events screen.
4. To change the date or time for the entry, tap **Time** and make the appropriate changes.
5. Use the number pad to enter the insulin amount.
6. Make sure the value displayed above the number pad is correct. If not, tap  to clear it and enter the correct value.
7. Tap **Log** at the bottom of the screen. The app returns to the Home screen, and  appears on the graph at the selected time.

Entering carbohydrate information

Use the app to record information about the carbohydrates consumed.

To enter carbohydrate information:

1. Determine the total amount, in grams, of carbohydrates in the planned meal, snack, or drink.
2. At the bottom of the Home screen, tap .
3. On the Events screen, tap .
4. To change the date or time for the entry, tap **Time** and make the appropriate changes.
5. Use the number pad to enter the amount of carbohydrates consumed.
6. Make sure the value displayed above the number pad is correct. If not, tap  to clear it and enter the correct value.
7. Tap **Log** at the bottom of the screen. The app returns to the Home screen, and  appears on the graph at the selected time.

Entering exercise information

Use the app to enter information about exercise regimen. Make sure to be consistent and enter the marker either before or after each period of exercise.

To enter exercise information:

1. Make a note about the length of time exercised (duration) and how difficult or easy the exercise was (intensity).
2. At the bottom of the Home screen, tap **+**.
3. On the Events screen, tap **+**.
4. To change the date, time, or duration, tap **Time** and set the desired day and time for this exercise entry. After tapping **OK**, set the exercise duration and tap **OK** to save the changes.
5. Tap **High**, **Medium**, or **Low** to indicate the intensity of exercise.
6. Tap **Log** at the bottom of the screen. The app returns to the Home screen, and **+** appears on the graph at the selected time.

Entering notes

Use the app to enter events other than BG measurements, insulin injections, carbohydrates consumed, and exercise information. For example, record information such as when medications are taken, illness, or stress.

To enter notes:

1. At the bottom of the Home screen, tap **+**.
2. On the Events screen, tap **+**.
3. Use the text field to enter the relevant information.
4. To change the date or time, tap **Time**. Tap **OK** to save the changes.
5. Tap **Save** at the top of the screen. The app returns to the Home screen, and **+** appears on the graph at the selected time.

Before using SG readings to make treatment decisions

Before using SG readings to make treatment decisions, consult with the healthcare professional to do the following:

- Develop a diabetes management plan
- Determine personal glucose target ranges

If the SG readings do not match symptoms, use a BG meter reading to confirm the SG value. If SG readings continue to be different than symptoms, consult the healthcare professional about how to use SG readings to help manage diabetes.

When to use BG meter readings

In the following conditions, use the BG meter readings to make treatment decisions.

- **The medication taken contains paracetamol or acetaminophen.**
Wait to use the SG reading to make treatment decisions until no longer taking medication. Medication that contains paracetamol or acetaminophen used to reduce fevers or treat cold and flu symptoms, can falsely raise SG readings.
- **The most recent SG reading is unavailable.**

If a new sensor is inserted or the Sensor updating message is displayed the SG reading will be unavailable. Check the BG value with a BG meter reading, and use the BG reading to make treatment decisions until sensor readings are available.

- **Symptoms do not match the SG value.**

Check the BG value with a BG meter reading before using the SG value to make treatment decisions.

CAUTION: Only use blood from fingertips to check BG levels. Do not use any other part of the body to obtain blood to test BG.

- Always make sure hands are clean before testing BG.
- Avoid use of an old BG reading or reuse of BG readings.

Note: The Simplera app contains a feature which provides alerts prior to SG levels reaching a set high or low limit. SG readings can vary from actual BG readings, therefore there may be situations when alerts are displayed and BG levels have not reached the set high or low limit. Consult the healthcare professional for actions and adjustments to alerts. There may also be situations when BG levels have reached the set high or low limit without any alerts displayed on the Simplera app. If there are symptoms of high or low BG levels, check the BG value immediately. Consult the healthcare professional for actions and adjustments for high or low BG levels.

Using SG readings to make treatment decisions

When using CGM, there are several things to consider to help make treatment decisions.

Look at the most recent SG value along with the SG graph, trend arrows, and alerts. The SG graph helps to understand how SG values may have recently changed. The trend arrows indicate what the SG level may be soon.

- If the SG reading is lower than the SG target, and CGM shows **↑↑**, this indicates the SG value is increasing. As a result, consider waiting to treat or adjust treatment for the low SG value. Consider symptoms before making treatment decisions based on the SG value.
- If the SG reading is higher than the SG target, and CGM shows **↓↓**, this indicates the SG value is decreasing. As a result, consider waiting to treat or adjust treatment for the high SG value. Consider symptoms before making treatment decisions based on the SG value.

Before taking an insulin dose based on the current SG, consider whether insulin from a previous dose is still lowering glucose levels. Stacking insulin is the process of taking an additional dose of insulin while there is still active insulin. Insulin stacking may cause low BG levels.

The following table can help to make treatment decisions.

Arrows Displayed	Low glucose	High glucose	Target glucose
None	May need to take a fast-acting sugar.	May adjust insulin to correct a high glucose to reach the glucose target range. Do not stack insulin.	Consider last insulin dose and carbs taken. Watch for any changes to glucose levels. Do not stack insulin.
1 UP Arrow ↑	Watch for any changes to glucose levels, con-	May adjust insulin to correct a high glucose to reach the glucose target range.	Consider last insulin dose and carbs taken. May need to take insulin

Arrows Displayed	Low glucose	High glucose	Target glucose
	sider symptoms, and wait to see if SG returns to target range. Make sure to not overtreat for a low SG level.	Do not stack insulin.	to stay within the glucose target range. Do not stack insulin.
2 UP Arrows ↑↑	Watch for any changes to glucose levels, consider symptoms, and wait to see if SG returns to target range. Make sure to not overtreat for a low SG level.	May adjust insulin to correct a high glucose to reach the glucose target range. Do not stack insulin.	May need to take insulin to stay within the glucose target range if insulin was not already taken with a recent meal or snack. Do not stack insulin.
3 UP arrows ↑↑↑	Watch for any changes to glucose levels, consider symptoms, and wait to see if SG returns to target range. Make sure to not overtreat for a low SG level.	May adjust insulin to correct a high glucose to reach the glucose target range. Do not stack insulin.	May need to take insulin to stay within target range. Do not stack insulin.
1 DOWN arrow ↓	May need to take a fast-acting sugar or eat a snack.	Consider last insulin dose and activity taken. May need to watch and wait to reach target range. Do not stack insulin.	May need to take a fast-acting sugar and eat a snack.
2 DOWN arrows ↓↓	May need to take a fast-acting sugar.	Consider last insulin dose and activity taken. May need to watch and wait to reach target range. Do not stack insulin.	May need to take a fast-acting sugar and eat a snack.
3 DOWN arrows ↓↓↓	May need to take a fast-acting sugar.	Consider last insulin dose and activity taken. May need to watch and wait to reach target range. Do not stack insulin.	May need to take a fast-acting sugar and eat a snack.

Syncing data to the CareLink Personal website

The Simplera system allows daily uploads to the CareLink Personal website. This is done using an automatic feature called Sync to CareLink.

This automatic Sync to CareLink feature sends sensor history information used to create CareLink Personal reports. This information can be viewed by a patient on the CareLink Personal website at carelink.minimed.eu. Please note that the mobile device must be connected to the Internet to

send data to the website. If using a cellular connection, the mobile device provider's data rates may apply.

If the **Sync to CareLink** toggle is turned off, the app no longer sends sensor information to the CareLink Personal website.

The **Upload Now** button gives the option to immediately send sensor history data to the website for generating reports.

Sharing CareLink Personal data with care partners

Now that the app is synced with the CareLink Personal website, users can share their data with a care partner, such as a family member or friend. This data is sent approximately every five minutes when an internet connection is available. Care partners can visit carelink.minimed.eu on a personal computer to create an account. Use the Manage Care Partners feature to give another person access to your information.

For more information, refer to the *CareLink Personal User Guide* that can be found on the CareLink Personal website.

Removing the Simplera sensor

To remove the Simplera sensor:

1. Gently peel the sensor adhesive away from the body.
2. Dispose of the Simplera sensor in accordance with all local laws and regulations. For additional information, see *Disposal, page 47*.

Bathing and swimming

While on the body, the Sensor is protected against continuous immersion in water at a depth of 2.4 meters (8 feet) for up to 30 minutes. Shower and swim without removing the sensor.

Troubleshooting

The following table contains troubleshooting information for the alerts.

Alerts

Problem	Likely Cause(s)	Resolution
Lost communication alert	Another app is in use, such as a game, that takes up a lot of the memory on the mobile device. This means that the Simplera app stops running and can't communicate with the sensor.	Open the app to ensure it is running properly. Check periodically to see if the app is still running in the background in order to receive alerts and SG values.
	The app has been closed. The Simplera app has stopped running and can't communicate with the sensor. Potential causes include using other apps and features (for example task manager apps), or selecting Force Stop for the Simplera app from an Android device Settings menu.	
Lost communication alert	The mobile device is out of range.	Make sure the mobile device and the sensor are located within 6 meters (20 feet). It is helpful to keep the devices on the same side of the body to minimize any radio frequency (RF) interference.
Lost communication alert	There is RF interference from other devices.	Move away from any equipment that can cause RF interference, such as cordless phones or routers.
Lost communication alert	The sensor pulled out from the skin.	Make sure that the sensor is still inserted. If the sensor has been pulled out from the skin, a new sensor must be used.
Mobile device battery low alert	The mobile device battery level is at 20% or lower. The battery needs to be recharged soon.	Recharge the mobile device battery to ensure that the Simplera app can function and send alerts. Remember to always carry a charger for the mobile device to ensure continuous use of the Simplera app.
Change sensor alert	The current sensor does not work properly and needs to be replaced.	To continue to receive SG values, a new sensor must be used. See <i>Inserting the Simplera sensor</i> , page 22 for instructions on how to insert the sensor.
Sensor end of life alert	The current sensor has reached the end of its life and will no longer display SG values on the Simplera app.	To continue to receive SG values, a new sensor must be used. See <i>Inserting the Simplera sensor</i> , page 22 for instructions on how to insert the sensor.

Problem	Likely Cause(s)	Resolution
Calibration not accepted alert	The last calibration value entered was not accepted by the Simplera sensor.	Wait at least 15 minutes before attempting another calibration. The Simplera app will request another calibration, if needed, after 15 minutes from when the Calibration not accepted alert was received. Review the BG meter instructions for use on how to test BG. Enter this new value in the app for calibration.
Connection error alert	There is an error with the sensor.	No action is required. The sensor is updating and this may take 30 minutes. During this time, do not rely on alerts from the app as SG information will not be available. Monitor BG levels with the meter.
Sensor updating alert	There is an error with the sensor.	No action is required. The sensor is updating and this may take up to 2 hours. During this time, do not rely on alerts from the app as SG information will not be available. Monitor BG levels with the meter.
Jailbroken device detected alert	The software on the iOS mobile device has been changed so that it no longer works in the way that the manufacturer intended.	The app cannot be used on a jailbroken device. To use the Simplera app the device software must remain as the manufacturer has designed. Changing the mobile device operating software causes the app to stop working.
Rooted device detected alert	The software on the Android mobile device has been changed so that it no longer works in the way that the manufacturer intended.	The app cannot be used on a rooted device. To use the Simplera app the device software must remain as the manufacturer has designed. Changing the mobile device operating software causes the app to stop working.

CareLink

Problem	Likely Cause(s)	Resolution
Unable to upload data to the CareLink Personal website	The phone has lost connection to the internet. The password has been changed through the CareLink website. CareLink servers are temporarily not responding.	<ol style="list-style-type: none"> 1. Check that the phone is connected to a Wi-Fi or mobile network data. Go to a web browser and try to access any web page to see if the internet connection is working through Wi-Fi or through mobile network data. 2. If the internet can be accessed, go to Menu > CareLink screen and re-log into CareLink again with latest CareLink credentials. 3. If credentials are updated, there is internet connection, and upload is still failing, CareLink servers may be temporarily not responding. Please try again later.

Specifications

Biocompatibility	Sensor: Complies with EN ISO 10993-1
Applied parts	Sensor
Operating conditions	Sensor temperature: 2 °C to 40 °C (36 °F to 104 °F) Sensor relative humidity: 15% to 95% with no condensation Sensor pressure: 70.33 kPa to 106.17 kPa (10.2 psi to 15.4 psi)
Storage conditions	Sensor temperature: 2 °C to 30 °C (36 °F to 86 °F) Sensor relative humidity: up to 95% with no condensation Sensor pressure: 70.33 kPa to 106.17 kPa (10.2 psi to 15.4 psi) CAUTION: Do not freeze the Simplera sensor, or store it in direct sunlight, extreme temperatures, or high humidity. These conditions may damage the device.
Duration of use	Up to seven days of CGM immediately following insertion.
Operating frequency	2.4 GHz band, Bluetooth® wireless technology (version 4.2)
Effective radiated power (ERP)	1.53 mW (1.85 dBm)
Effective isotropic radiated power (EIRP)	2.51 mW (4.00 dBm)
Operating range	Up to 6.09 m (20 ft) line of sight in free-air
Sensor device approximate dimensions	6.009 x 6.009 x 7.414 cm (2.366 x 2.366 x 2.919 in)
Sensor device approximate weight	72.5 g (2.56 ounces)
Sensor approximate dimensions	2.865 x 2.865 x 0.477 cm (1.128 x 1.128 x 0.188 in)
Sensor approximate weight	4.6 g (0.16 ounces)

Maintenance

Cleaning

Not applicable.

Disinfecting

Not applicable.

Disposal

Disposal requirements for electronic equipment, batteries, sharps and potential biohazardous materials differ based on location. Confirm disposal requirements for electronic equipment, batteries, sharps, and potential biohazardous materials with local laws and regulations.

The used inserter contains a needle which has been in contact with blood or other bodily fluids.

The used sensor contains a battery and has been in contact with blood or other bodily fluids. Disposal of the battery in any receptacle that could be exposed to extreme heat may cause the battery to ignite and result in serious injury.

Do not dispose of any component of this product with household waste or recyclables.

Dispose of the inserter and sensor in accordance with local laws and regulations.

Life of use

The Simplera sensor can be used one time and has a maximum life of up to 170 hours (seven days). The 170-hour life span of the sensor begins within 30 minutes after insertion.

CAUTION: Do not use the sensor if there is a sudden rise in sensor temperature. When operating the sensor in air temperatures of 40 °C (104 °F), under certain fault conditions, the temperature of the sensor may briefly rise up to 50 °C (121 °F). If there is a sudden rise in temperature or the sensor becomes hot or uncomfortable, remove and discard the sensor.

Simplera app quality of service

The Simplera app can use either Wi-Fi or cellular data to send data to the CareLink Connect app for remote monitoring, and to upload history to the CareLink Personal website. The Simplera app will use Wi-Fi to transmit data when a Wi-Fi connection is available, and cellular data if Wi-Fi is not available. Although all data sent by the Simplera app is encrypted, a secured Wi-Fi network is recommended.

The Simplera sensor connects to a compatible mobile device via a Bluetooth low energy technology network. The sensor sends glucose data and status-related alerts to the compatible mobile device, which verifies the integrity of received data after wireless transmission. The quality of the connection is in accordance with the Bluetooth Specification v4.2.

Data security

The Simplera sensor is designed to only accept radio frequency (RF) communications from a recognized and linked compatible mobile device. The sensor must be paired with the mobile device before the mobile device accepts information from the sensor.

The compatible mobile device ensures data security via proprietary means and data integrity using error checking processes, such as cyclic redundancy checks.

The Simplera app has been designed with security features to help keep its data secure. However, there are important recommended steps to take to ensure the compatible mobile device used with the Simplera app is also secure:

- Do not leave the compatible mobile device unattended.
- Use caution when viewing or sharing data with others.
- Enable a security lock on the compatible mobile device. When the compatible mobile device is not in use, lock it in a way that requires the password to be entered in order to use it.
- Do not remove or interfere with the security features on the compatible mobile device.
- Do not attempt to modify the operating system, jailbreak, or root the device.
- Use only the official application store, such as the Apple App Store or the Google Play store to get all mobile applications used with the compatible mobile device.
- Do not click on links from email messages, web pages, or text messages received from an unknown or untrusted source.
- Avoid the use of unknown Wi-Fi networks or public Wi-Fi hotspots.
- Enable security protection on a home Wi-Fi network, such as the use of a password and encryption.

The app may send anonymous analytic data to Medtronic if permission has been granted in the setup of the app. This data is used to analyze crash logs and app performance. This access can be revoked or reinstated at any time in the Carelink screen of the app.

Traveling by air

The Simplera sensor is safe for use on commercial airlines.

Open Source Software (OSS) disclosure

This document identifies the Open Source Software that may be separately called, executed, linked, affiliated, or otherwise utilized by this product.

Such Open Source Software is licensed to users subject to the terms and conditions of the separate software license agreement for such Open Source Software.

Use of the Open Source Software by you shall be governed entirely by the terms and conditions of such license.

The source/object code and applicable license for the Open Source Software can be obtained at the following site: www.medtronicdiabetes.com/ossnotices.

Guidance and manufacturer's declaration

Guidance and Manufacturer's Declaration - Electromagnetic Emissions		
Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions CISPR 11	CISPR 11 Group 1, Class B	The transmitter uses RF energy only for system communications. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
Immunity Test	IEC 60601-1-2 Test Level	IEC 60601-1-2 Compliance Level	Electromagnetic Environment Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	For use in a typical domestic, commercial, or hospital environment.
Power frequency magnetic field IEC 61000-4-8	30 A/m	30 A/m	For use in a typical domestic, commercial, or hospital environment.
Proximity magnetic fields IEC 61000-4-39, Table 11	IEC 60601-1-2, Table 11	IEC 60601-1-2, Table 11	For use in a typical domestic, commercial, or hospital environment.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity

Immunity Test	IEC 60601-1-2 Test Level	IEC 60601-1-2 Compliance Level	Electromagnetic Environment Guidance
Proximity fields from RF wireless communications equipment	IEC 60601-1-2, Table 9	IEC 60601-1-2, Table 9	For use in a typical domestic, commercial, or hospital environment. Portable and mobile RF communications equipment should be used no closer to any part of the transmitter than the recommended separation distance of 30 cm (12 in).
Radiated RF electromagnetic fields IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz 80% AM at 1 kHz	10 V/m 80 MHz to 2.7 GHz 80% AM at 1 kHz	Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol: 
<p>Note: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption, and reflection from structures, objects and people.</p>			

Simplera system performance

Note: A patient should review the information in this section with a healthcare professional to understand the performance of the Simplera system.

Clinical study overview

The performance of the Simplera system was evaluated using data collected during a multi-center prospective clinical study.¹ The study included participants ages 2 to 80 years old. The study enrolled a total of 243 subjects previously diagnosed with type 1 or 2 diabetes and 230 of these subjects completed the study. Subjects 18 and older were instructed to wear a total of two sensors, in the arm. Subjects age 7 to 17 were instructed to wear a total of three sensors in the arm(s) and buttock(s). Subjects age 2 to 6 were instructed to wear a total of two sensors in the arm and buttock. For all subjects, the sensors were used to record raw sensor signals during the study and there was no real-time calculation of sensor glucose values.

¹ NCT04436822: CIP330 - Evaluation of Updated Continuous Glucose Monitoring (CGM) Form Factor in Adults, Adolescents and Pediatrics.

Frequent sample testing (FST) was performed on four occasions for subjects 14 and older and on two occasions for subjects 2 to 13 years of age.

During FST, reference blood (plasma) glucose values were obtained with a Yellow Springs Instrument (YSI™) Glucose Analyzer every 5-15 minutes for subjects 7 and older. Fingerstick blood glucose meter measurements were used as reference values for pediatric subjects aged 2 to 6 years old. During each FST, subjects 14 and older with an established insulin sensitivity ratio and insulin carbohydrate ratio underwent a hypoglycemic challenge or a hyperglycemic challenge.

Data collected during the study was post-processed after the study using the Simplera system sensor algorithm to convert the raw sensor information to sensor glucose values every five minutes. For the accuracy information presented in the following section, YSI or finger stick reference values were paired with the closest sensor glucose reading within five minutes of the time of the reference value measurement.

Sensor accuracy

Sensor accuracy was calculated for sensors compared to a self-monitoring blood glucose (SMBG) reference for 2-6 year old subjects and compared to a YSI reference for subjects ages 7 and older.

Table 1. Overall Accuracy Compared to YSI

Patient Population	Insertion Site	Number of Subjects	Number of paired SG-YSI	Percent of SG within 20%** of YSI (95% Lower Bound)**	Mean Absolute Relative Difference (%)
Adults (18+)	Arm	116	15405	90.7 (90.3)	10.2
	Arm	111	8632	88.6 (88.0)	10.9
Pediatrics (2-17)*	Buttocks	99	7783	89.6 (89.0)	10.2
	Arm	89	8282	89.0 (88.4)	10.8
Pediatrics (7-17)	Buttocks	86	7569	89.9 (89.3)	10.1
	Arm	22	350	79.7 (76.2)	13.2
Pediatrics (2-6)*	Buttock	13	214	79.0 (74.4)	13.6

*SMBG were collected for subjects ages 2-6 years and were used in place of YSI measurements.
**For 20% agreement, 20 mg/dL used when YSI <70 mg/dL.

CGM readings are within 50-400 mg/dL, inclusive.

Table 2. Simplera system accuracy, summary of Mean Absolute Relative Difference (MARD) by FST day, Adults, Arm insertion site

Reference Range	Characteristic	FST1	FST2	FST3	FST4	Total
<3.9 mmol/L (<70 mg/dL)	No. of Paired Points	786	812	448	219	2265
	Mean (SD)	17.7 (14.5)	12.6 (11.9)	9.6 (8.5)	12.6 (12.6)	13.7 (12.8)
	Median	15.4	8.8	7.8	9.2	10.4
3.9-10 mmol/L (70-180 mg/dL)	No. of Paired Points	2488	2586	2077	1667	8818
	Mean (SD)	11.0 (10.7)	9.9 (8.1)	9.5 (7.5)	10.0 (8.5)	10.1 (8.9)
	Median	7.9	8	7.9	8.2	8
>10 mmol/L (>180 mg/dL)	No. of Paired Points	1096	846	1345	1035	4322
	Mean (SD)	10.5 (8.0)	8.3 (6.0)	7.8 (6.5)	8.2 (6.7)	8.7 (7.0)
	Median	9.3	7.1	6.5	6.8	7.3
Overall	No. of Paired Points	4370	4244	3870	2921	15405
	Mean (SD)	12.0 (11.2)	10.1 (8.7)	8.9 (7.3)	9.5 (8.4)	10.2 (9.2)
	Median	9	7.9	7.4	7.7	8

Table 3. Simplera system accuracy, summary of MARD by FST day, Pediatrics, (7-17), Arm insertion site

Reference Range	Characteristic	FST1	FST2	FST3	FST4	Total
<3.9 mmol/L (<70 mg/dL)	No. of Paired Points	479	391	270	240	1380
	Mean (SD)	21.3 (15.3)	12.7 (17.5)	12.3 (9.2)	9.5 (8.4)	15.1 (14.8)
	Median	18.7	7.4	10.2	6.5	11.1
3.9-10 mmol/L (70-180 mg/dL)	No. of Paired Points	1478	1345	584	467	3874
	Mean (SD)	11.5 (9.9)	10.8 (10.4)	9.4 (8.6)	9.1 (9.1)	10.6 (9.8)
	Median	9.1	8.1	6.8	6.6	8.1
>10 mmol/L (>180 mg/dL)	No. of Paired Points	952	908	759	409	3028
	Mean (SD)	9.8 (7.8)	9.8 (9.7)	7.2 (6.1)	9.2 (8.3)	9.1 (8.2)
	Median	7.6	7.4	5.6	7.3	7

Table 3. Simpler system accuracy, summary of MARD by FST day, Pediatrics, (7-17), Arm insertion site (continued)

Reference Range	Characteristic	FST1	FST2	FST3	FST4	Total
Overall	No. of Paired Points	2909	2644	1613	1116	8282
	Mean (SD)	12.5 (11.2)	10.7 (11.6)	8.9 (7.8)	9.2 (8.7)	10.8 (10.5)
	Median	9.8	7.8	6.8	6.8	8

Includes pediatric subjects 7-17 years of age; no YSI measurements were collected for pediatric subjects 2-6 years of age.

Table 4. Simpler system accuracy, summary of MARD by FST day, Pediatrics, (7-17), Buttock insertion site

Reference Range	Characteristic	FST1	FST2	FST3	FST4	Total
<3.9 mmol/L (<70 mg/dL)	No. of Paired Points	410	371	272	202	1255
	Mean (SD)	22.2 (13.2)	16.2 (18.9)	9.2 (8.6)	13.2 (9.5)	16.2 (14.7)
	Median	21.1	10.6	6.6	13.2	12.7
3.9-10 mmol/L (70-180 mg/dL)	No. of Paired Points	1368	1260	537	342	3507
	Mean (SD)	11.3 (10.0)	9.6 (9.7)	8.1 (7.5)	8.4 (6.7)	10.0 (9.3)
	Median	8.7	7	6.4	7	7.6
>10 mmol/L (>180 mg/dL)	No. of Paired Points	932	909	642	324	2807
	Mean (SD)	7.9 (6.8)	7.3 (6.9)	8.2 (7.6)	6.9 (6.3)	7.6 (7.0)
	Median	6.1	5.2	6.2	5.2	5.7
Overall	No. of Paired Points	2710	2540	1451	868	7569
	Mean (SD)	11.8 (10.7)	9.7 (11.1)	8.3 (7.8)	9.0 (7.7)	10.1 (10.1)
	Median	8.7	6.7	6.4	6.8	7.3

Includes pediatric subjects 7-17 years of age; no YSI measurements were collected for pediatric subjects 2-6 years of age.

Alert performance

Glucose true alert rate

The glucose true alert rate is the rate at which the BG confirmed that the CGM alert was triggered correctly.

Table 5. Glucose True Alert Rate, Adults

Glucose	Insertion Site	Glucose TRUE Alert Rate					
		Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
2.8 mmol/L (50 mg/dL)	Arm	31.8%	31.8%	14.9%	13.3%	17.5%	16.1%
3 mmol/L (54 mg/dL)	Arm	51.7%	51.7%	24.5%	20.4%	29.2%	25.9%
3.3 mmol/L (60 mg/dL)	Arm	72.2%	71.3%	36.4%	32.0%	45.5%	42.0%
3.9 mmol/L (70 mg/dL)	Arm	81.0%	80.4%	42.5%	41.0%	53.3%	52.0%
4.4 mmol/L (80 mg/dL)	Arm	79.4%	78.0%	44.9%	41.8%	55.4%	52.8%
5 mmol/L (90 mg/dL)	Arm	75.9%	75.9%	49.2%	46.1%	58.5%	56.4%
10 mmol/L (180 mg/dL)	Arm	88.5%	88.3%	63.0%	60.5%	72.5%	70.8%
12.2 mmol/L (220 mg/dL)	Arm	89.8%	89.1%	60.8%	58.5%	71.2%	69.5%
13.9 mmol/L (250 mg/dL)	Arm	90.1%	89.5%	57.7%	55.2%	68.5%	66.7%
16.7 mmol/L (300 mg/dL)	Arm	95.7%	95.7%	62.0%	57.2%	72.5%	69.2%

Table 6. Glucose True Alert Rate, Pediatrics, (7-17)

Glucose (mmol/L)	Insertion Site	Glucose True Alert Rate					
		Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
2.8	Buttock	48.1%	48.1%	15.3%	11.7%	20.0%	16.8%
	Arm	52.4%	52.4%	15.1%	13.7%	19.8%	18.6%

Table 6. Glucose True Alert Rate, Pediatrics, (7-17) (continued)

Glucose True Alert Rate							
Glucose (mmol/L)	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
3	Buttock	71.0%	71.0%	20.3%	16.5%	27.7%	24.4%
	Arm	65.6%	65.6%	19.9%	17.6%	26.9%	25.0%
3.3	Buttock	85.7%	85.7%	37.1%	32.7%	46.6%	43.0%
	Arm	69.1%	69.1%	33.3%	27.4%	41.4%	36.8%
3.9	Buttock	78.7%	78.7%	45.5%	41.5%	54.7%	51.8%
	Arm	75.7%	74.8%	42.4%	38.5%	51.4%	48.4%
4.4	Buttock	78.6%	77.8%	47.7%	44.7%	57.7%	55.4%
	Arm	71.9%	71.2%	45.2%	43.0%	54.1%	52.4%
5	Buttock	83.3%	82.1%	52.1%	49.0%	63.2%	60.8%
	Arm	76.5%	76.0%	50.9%	48.2%	59.8%	57.8%
10	Buttock	88.8%	88.8%	68.9%	66.3%	76.8%	75.3%
	Arm	89.5%	89.1%	73.2%	69.8%	79.9%	77.8%
12.2	Buttock	93.2%	93.2%	67.3%	64.0%	77.0%	75.0%
	Arm	93.9%	93.4%	68.6%	65.4%	78.2%	76.0%
13.9	Buttock	91.3%	90.6%	61.3%	56.4%	72.2%	68.8%
	Arm	90.3%	89.7%	63.8%	59.0%	73.4%	70.1%
16.7	Buttock	88.1%	88.1%	52.7%	50.0%	65.2%	63.4%
	Arm	89.7%	89.7%	57.2%	55.0%	68.6%	67.1%

Includes pediatric subjects 7-17 years of age; no YSI measurements were collected for pediatric subjects 2-6 years of age.

Glucose false alert rate

The glucose false alert rate is the rate at which the BG did not confirm that the CGM alert was triggered correctly.

Table 7. Glucose False Alert Rate, Adults

Glucose FALSE Alert Rate							
Glucose	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
2.8 mmol/L (50 mg/dL)	Arm	68.2%	68.2%	85.1%	86.7%	82.5%	83.9%
3 mmol/L (54 mg/dL)	Arm	48.3%	48.3%	75.5%	79.6%	70.8%	74.1%
3.3 mmol/L (60 mg/dL)	Arm	27.8%	28.7%	63.6%	68.0%	54.5%	58.0%
3.9 mmol/L (70 mg/dL)	Arm	19.0%	19.6%	57.5%	59.0%	46.7%	48.0%
4.4 mmol/L (80 mg/dL)	Arm	20.6%	22.0%	55.1%	58.2%	44.6%	47.2%
5 mmol/L (90 mg/dL)	Arm	24.1%	24.1%	50.8%	53.9%	41.5%	43.6%
10 mmol/L (180 mg/dL)	Arm	11.5%	11.8%	37.0%	39.5%	27.5%	29.2%
12.2 mmol/L (220 mg/dL)	Arm	10.2%	10.9%	39.2%	41.5%	28.8%	30.5%
13.9 mmol/L (250 mg/dL)	Arm	9.9%	10.5%	42.3%	44.8%	31.5%	33.3%
16.7 mmol/L (300 mg/dL)	Arm	4.3%	4.3%	38.0%	42.8%	27.5%	30.8%

Table 8. Glucose False Alert Rate, Pediatrics, (7-17)

Glucose False Alert Rate							
Glucose (mmol/L)	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
2.8	Buttock	51.9%	51.9%	84.7%	88.3%	80.0%	83.2%
	Arm	47.6%	47.6%	84.9%	86.3%	80.2%	81.4%
3	Buttock	29.0%	29.0%	79.7%	83.5%	72.3%	75.6%
	Arm	34.4%	34.4%	80.1%	82.4%	73.1%	75.0%
3.3	Buttock	14.3%	14.3%	62.9%	67.3%	53.4%	57.0%
	Arm	30.9%	30.9%	66.7%	72.6%	58.6%	63.2%
3.9	Buttock	21.3%	21.3%	54.5%	58.5%	45.3%	48.2%
	Arm	24.3%	25.2%	57.6%	61.5%	48.6%	51.6%
4.4	Buttock	21.4%	22.2%	52.3%	55.3%	42.3%	44.6%
	Arm	28.1%	28.8%	54.8%	57.0%	45.9%	47.6%

Table 8. Glucose False Alert Rate, Pediatrics, (7-17) (continued)

Glucose False Alert Rate							
Glucose (mmol/L)	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
5	Buttock	16.7%	17.9%	47.9%	51.0%	36.8%	39.2%
	Arm	23.5%	24.0%	49.1%	51.8%	40.2%	42.2%
10	Buttock	11.2%	11.2%	31.1%	33.7%	23.2%	24.7%
	Arm	10.5%	10.9%	26.8%	30.2%	20.1%	22.2%
12.2	Buttock	6.8%	6.8%	32.7%	36.0%	23.0%	25.0%
	Arm	6.1%	6.6%	31.4%	34.6%	21.8%	24.0%
13.9	Buttock	8.8%	9.4%	38.7%	43.6%	27.8%	31.2%
	Arm	9.7%	10.3%	36.2%	41.0%	26.6%	29.9%
16.7	Buttock	11.9%	11.9%	47.3%	50.0%	34.8%	36.6%
	Arm	10.3%	10.3%	42.8%	45.0%	31.4%	32.9%

Includes pediatric subjects 7-17 years of age; no YSI measurements were collected for pediatric subjects 2-6 years of age.

Glucose correct detection alert rate

Glucose Correct Detection Rate is the rate that the device alerted when it should have alerted.

Table 9. Glucose Correct Detection Alert Rate, Adults

Glucose Correct Detection Alert Rate							
Glucose	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
2.8 mmol/L (50 mg/dL)	Arm	41.2%	41.2%	97.1%	91.2%	97.1%	91.2%
3 mmol/L (54 mg/dL)	Arm	47.6%	47.6%	96.8%	90.5%	96.8%	90.5%
3.3 mmol/L (60 mg/dL)	Arm	65.5%	64.7%	92.2%	84.5%	92.2%	84.5%
3.9 mmol/L (70 mg/dL)	Arm	82.2%	82.2%	94.7%	92.1%	94.7%	92.8%
4.4 mmol/L (80 mg/dL)	Arm	86.8%	84.9%	96.2%	93.4%	96.2%	93.9%
5 mmol/L (90 mg/dL)	Arm	84.9%	84.6%	94.4%	90.2%	94.4%	91.9%
10 mmol/L (180 mg/dL)	Arm	90.2%	88.3%	98.4%	95.1%	98.6%	95.8%
12.2 mmol/L (220 mg/dL)	Arm	85.9%	85.6%	95.8%	93.3%	96.1%	94.0%
13.9 mmol/L (250 mg/dL)	Arm	80.3%	80.3%	93.6%	90.6%	93.6%	91.6%
16.7 mmol/L (300 mg/dL)	Arm	74.6%	73.7%	94.9%	89.0%	94.9%	89.8%

Table 10. Glucose Correct Detection Alert Rate, Pediatrics, (7-17)

Glucose Correct Detection Alert Rate							
Glucose (mmol/L)	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
2.8	Buttock	52.0%	52.0%	100.0%	76.0%	100.0%	80.0%
	Arm	42.3%	42.3%	88.5%	80.8%	88.5%	80.8%
3	Buttock	53.5%	53.5%	90.7%	74.4%	90.7%	76.7%
	Arm	51.2%	51.2%	86.0%	74.4%	86.0%	74.4%
3.3	Buttock	57.5%	57.5%	92.5%	87.5%	92.5%	87.5%
	Arm	61.7%	59.3%	92.6%	81.5%	92.6%	82.7%
3.9	Buttock	80.4%	79.3%	98.9%	96.7%	98.9%	96.7%
	Arm	80.0%	79.0%	99.0%	96.0%	99.0%	96.0%
4.4	Buttock	91.4%	90.5%	100.0%	99.1%	100.0%	99.1%
	Arm	92.0%	91.2%	100.0%	99.2%	100.0%	100.0%
5	Buttock	91.0%	91.0%	97.2%	94.4%	97.2%	95.1%
	Arm	87.4%	86.8%	98.7%	95.6%	98.7%	96.2%
10	Buttock	94.4%	93.4%	98.3%	95.5%	98.3%	97.9%
	Arm	93.7%	91.7%	96.7%	93.4%	97.3%	96.7%
12.2	Buttock	92.9%	90.7%	98.7%	96.4%	98.7%	97.8%
	Arm	90.7%	89.8%	96.6%	94.9%	96.6%	95.3%
13.9	Buttock	88.3%	88.3%	95.3%	93.0%	95.9%	94.2%
	Arm	86.9%	85.8%	95.6%	94.0%	95.6%	95.1%

Table 10. Glucose Correct Detection Alert Rate, Pediatrics, (7-17) (continued)

Glucose Correct Detection Alert Rate							
Glucose (mmol/L)	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
16.7	Buttock	91.8%	91.8%	95.9%	93.8%	95.9%	94.8%
	Arm	82.1%	81.3%	89.3%	84.8%	90.2%	88.4%

Includes pediatric subjects 7-17 years of age; no YSI measurements were collected for pediatric subjects 2-6 years of age.

Glucose missed detection alert rate

Glucose Missed Detection Rate is the rate that the device did not alert when it should have alerted.

Table 11. Glucose Missed Detection Alert Rate, Adults

Glucose Missed Detection Alert Rate							
Glucose	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
2.8 mmol/L (50 mg/dL)	Arm	58.8%	58.8%	2.9%	8.8%	2.9%	8.8%
3 mmol/L (54 mg/dL)	Arm	52.4%	52.4%	3.2%	9.5%	3.2%	9.5%
3.3 mmol/L (60 mg/dL)	Arm	34.5%	35.3%	7.8%	15.5%	7.8%	15.5%
3.9 mmol/L (70 mg/dL)	Arm	17.8%	17.8%	5.3%	7.9%	5.3%	7.2%
4.4 mmol/L (80 mg/dL)	Arm	13.2%	15.1%	3.8%	6.6%	3.8%	6.1%
5 mmol/L (90 mg/dL)	Arm	15.1%	15.4%	5.6%	9.8%	5.6%	8.1%
10 mmol/L (180 mg/dL)	Arm	9.8%	11.7%	1.6%	4.9%	1.4%	4.2%
12.2 mmol/L (220 mg/dL)	Arm	14.1%	14.4%	4.2%	6.7%	3.9%	6.0%
13.9 mmol/L (250 mg/dL)	Arm	19.7%	19.7%	6.4%	9.4%	6.4%	8.4%
16.7 mmol/L (300 mg/dL)	Arm	25.4%	26.3%	5.1%	11.0%	5.1%	10.2%

Table 12. Glucose Missed Detection Alert Rate, Pediatrics, (7-17)

Glucose Missed Detection Alert Rate							
Glucose (mmol/L)	Insertion Site	Threshold Only		Predictive Only		Threshold & Predictive	
		30 min	15 min	30 min	15 min	30 min	15 min
2.8	Buttock	48.0%	48.0%	0.0%	24.0%	0.0%	20.0%
	Arm	57.7%	57.7%	11.5%	19.2%	11.5%	19.2%
3	Buttock	46.5%	46.5%	9.3%	25.6%	9.3%	23.3%
	Arm	48.8%	48.8%	14.0%	25.6%	14.0%	25.6%
3.3	Buttock	42.5%	42.5%	7.5%	12.5%	7.5%	12.5%
	Arm	38.3%	40.7%	7.4%	18.5%	7.4%	17.3%
3.9	Buttock	19.6%	20.7%	1.1%	3.3%	1.1%	3.3%
	Arm	20.0%	21.0%	1.0%	4.0%	1.0%	4.0%
4.4	Buttock	8.6%	9.5%	0.0%	0.9%	0.0%	0.9%
	Arm	8.0%	8.8%	0.0%	0.8%	0.0%	0.0%
5	Buttock	9.0%	9.0%	2.8%	5.6%	2.8%	4.9%
	Arm	12.6%	13.2%	1.3%	4.4%	1.3%	3.8%
10	Buttock	5.6%	6.6%	1.7%	4.5%	1.7%	2.1%
	Arm	6.3%	8.3%	3.3%	6.6%	2.7%	3.3%
12.2	Buttock	7.1%	9.3%	1.3%	3.6%	1.3%	2.2%
	Arm	9.3%	10.2%	3.4%	5.1%	3.4%	4.7%
13.9	Buttock	11.7%	11.7%	4.7%	7.0%	4.1%	5.8%
	Arm	13.1%	14.2%	4.4%	6.0%	4.4%	4.9%
16.7	Buttock	8.2%	8.2%	4.1%	6.2%	4.1%	5.2%
	Arm	17.9%	18.8%	10.7%	15.2%	9.8%	11.6%

Includes pediatric subjects 7-17 years of age; no YSI measurements were collected for pediatric subjects 2-6 years of age.

Sensor life

Adults

96.6% of sensors worn in the arm operated more than six days and up to the full seven days of wear (144 to 168 hours). The mean functional sensor life for sensors worn in the arm insertion site over the course of the study was 147.8 hours, with a median functional life of 167.9 hours.

Pediatrics 7-17

94.7% of sensors worn in the arm operated more than six days and up to the full seven days of wear (144 to 168 hours). The mean functional sensor life for sensors worn in the arm insertion site over the course of the study was 146.3 hours, with a median functional life of 162.5 hours. 96.8% of sensors worn in the buttock operated more than six days and up to the full seven days of wear (144 to 168 hours). The mean functional sensor life for sensors worn in the buttock insertion site over the course of the study was 142.7 hours, with a median functional life of 161 hours.

Pediatrics 2-6

91.3% of sensors worn in the arm operated more than six days and up to the full seven days of wear (144 to 168 hours). The mean functional sensor life for sensors worn in the arm insertion site over the course of the study was 141.5 hours, with a median functional life of 167.9 hours. 76.9% of sensors worn in the buttock operated more than six days and up to the full seven days of wear (144 to 168 hours). The mean functional sensor life for sensors worn in the buttock insertion site over the course of the study was 151.2 hours, with a median functional life of 167.9 hours.

Safety

Device related adverse events were limited to pain or bruising at the sensor insertion site.

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© 2023 Medtronic
Medtronic B.V.
Medtronic Limited